Connecting BiWaze® Cough to Arc Connect™

Use this guide to set up the Arc Connect app and connect your BiWaze Cough device.

Before you start, have these ready:

- Your phone
- Wi-Fi name and password
- BiWaze Cough device turned on

1 Get the app

- Go to Apple App Store or Google Play
- Search for Arc Connect
- Download the app





Download on the App Store

GET IT ON Google Play

2 Create your account

- Open the app
- Tap “Register”
- Complete account creation steps

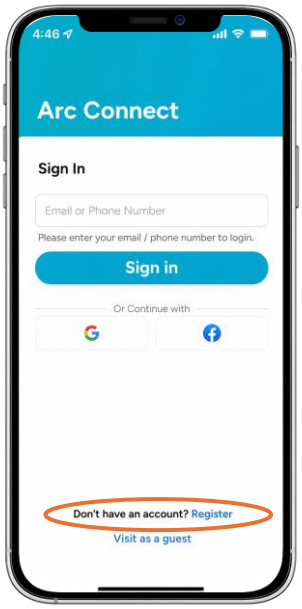


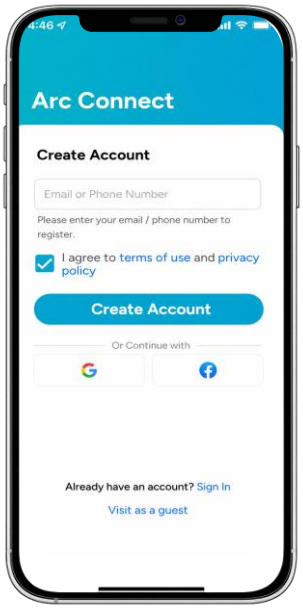

3 Connect to Wi-Fi on BiWaze Cough

- Open Wi-Fi
- Turn Wi-Fi on
- Choose network
- Enter Wi-Fi password



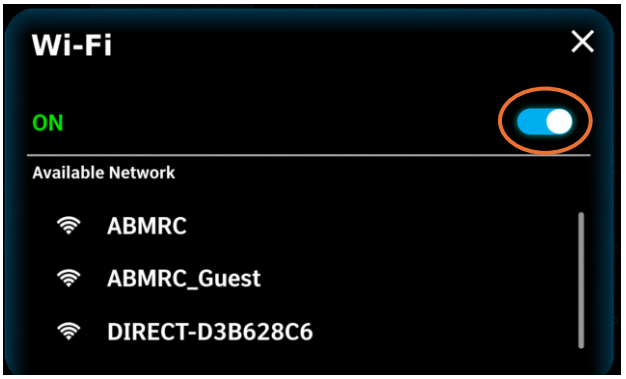


See page 2 to continue setup →

Connecting BiWaze® Cough to Arc Connect™

(Continued – page 2 of 2)

4 Open the QR code on BiWaze Cough


- Tap the gear icon on bottom left of home screen
- Make sure Arc Connect is ON
- Tap Arc Connect button
- Keep QR code on screen

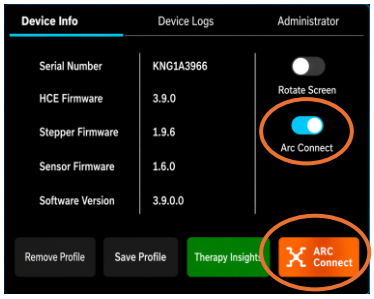




5 Link your device in the app

- Follow setup screens in app
- Select “BiWaze Cough”
- Tap “Continue”
- Tap “Scan the QR Code”

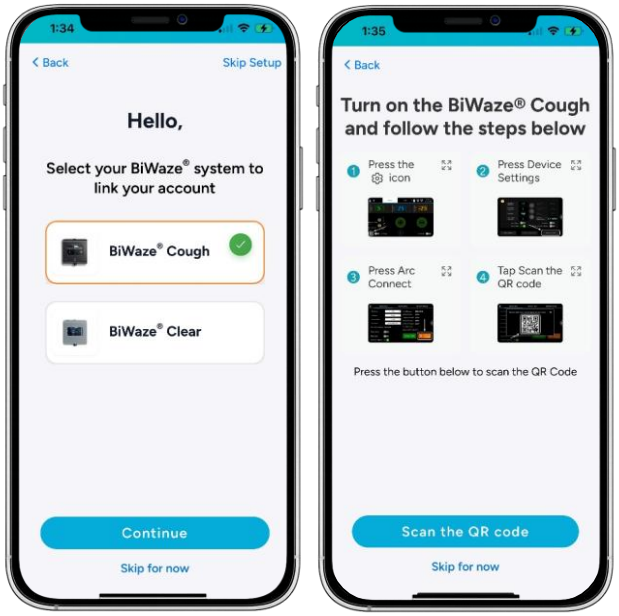

6 Finish setup

- Confirm device connection
- Set your goals and reminders, sharing preferences, and health app connection in My Settings (optional)

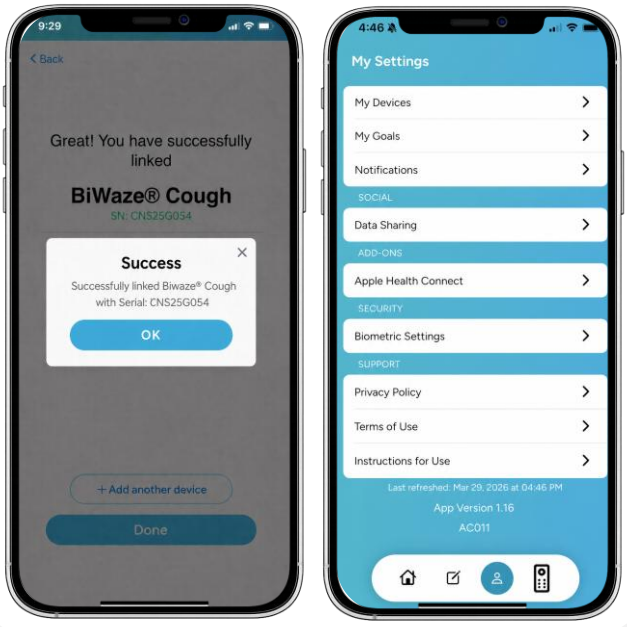

Need help?
Contact your Home Health team _____

BiWaze® and ABM Respiratory Care® are registered trademarks of ABMRC, LLC.
©2026 ABMRC, LLC. ALL RIGHTS RESERVED. PRTN-158 APR. 205539264-1492 Rev 1.0 MAY 2026
For further information about this product or a service, visit our webpage: www.mec.com.au or call customer support: 1300 632 633